


for MH6000P

Operation for return a factory setting:


1. Remove battery and AC/DC Power Supply and make sure remove **ALL** the electrodes from skin.
2. Waiting one minute.
3. To insert the battery into the compartment and make sure that you are installing the battery properly.

(make sure remove ALL the electrodes from skin)

4. **First** press and hold P+ Program number Increase Button and P- Program number Decrease buttons at the same time,

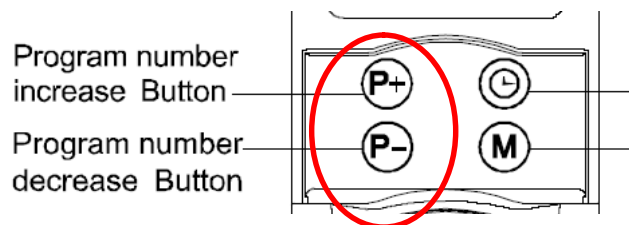

then Turn Intensity Control knob for Channel 1 or 2 clockwise for power on.

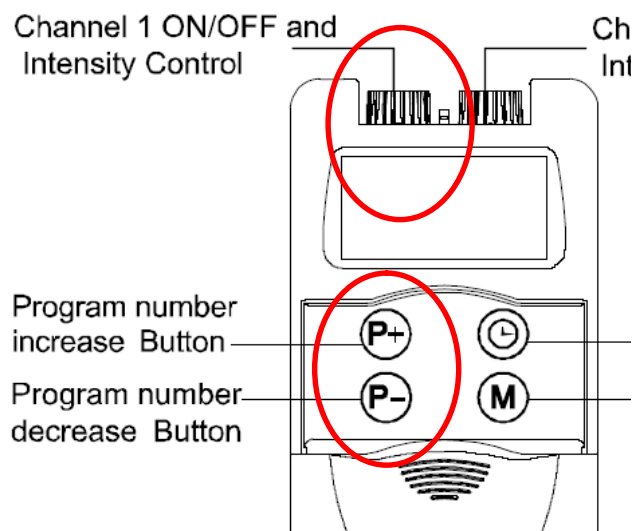

5. The unit will be signaling "Bi" sound about 2 seconds.
6. Please release all buttons.
7. Turn both Channel Intensity Controls to switch off unit.
8. The unit already returned a factory setting.