



Operation for return a factory setting:

1. Remove battery and make sure remove **ALL** the electrodes from skin.
2. Waiting one minute.
3. To insert the battery into the compartment and make sure that you are installing the battery properly.
(make sure remove **ALL** the electrodes from skin)
4. **First** press and hold two intensity decrease buttons,

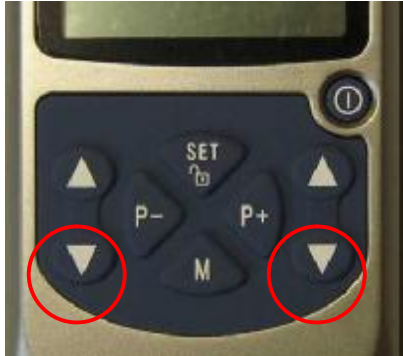

and press and hold power ON/OFF button for 1 second for power on.

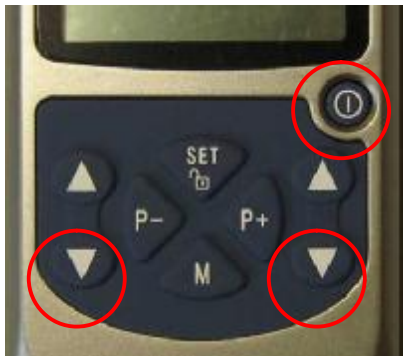

5. The unit will be signaling "Bi" sound about 2 seconds.
6. Please release all buttons.
7. Press the Power ON/OFF button to switch off unit.
8. The unit already returned a factory setting.